

OptiBond Solo Plus. Total-etch. Dual-cure. Great bonds.



BONDING AGENT

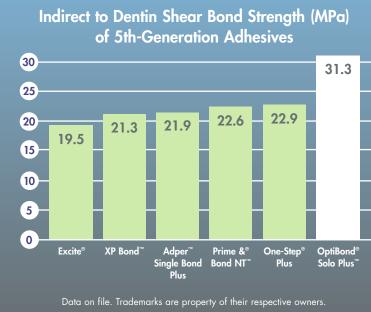

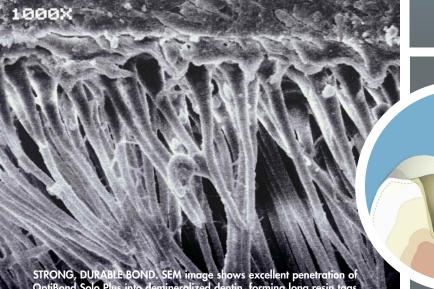

FEATURES

Filled technology Fluoride release Ethanol based Unidose delivery

BENEFITS

Strong bond Less polymerization shrinkage Less microleakage Less sensitivity Use in moist or dry environment

Four Star Award

2003-2006

The strength behind your restorations.

- Clinically proven. Universally accepted, industry-standard products. Proven performance worldwide.
- Strong, durable bond. 15%-filled material penetrates and reinforces dentin tubules to provide a long-lasting bond.
- **Delivery options.** Choose 5 ml bottle or convenient Unidose® delivery.
- Dual-cure option. OptiBond Dual Cure Activator provides added assurance of a total cure if the area is not accessible to adequate light.



OptiBond Solo Plus. Make our bond strength yours.

When it comes to choosing a clinically proven bonding agent, OptiBond Solo Plus is the gold standard, employing a single-component primer/adhesive in combination with Kerr's phosphoric acid gel etchant.

OptiBond Solo Plus is 15%-filled with the same 0.4 micron filler found in Kerr's industry-recognized composites. The filler reinforces the hybrid zone and the dentin tubules, creating a true structured bond of extraordinary strength not found in conventional, unfilled adhesive systems. OptiBond Solo Plus's durable chemical and micro-mechanical bonds protect against microleakage and decrease polymerization shrinkage to ensure superior marginal integrity.

Indications: Direct composite to enamel or dentin, composite to composite, composite to porcelain or metal, amalgam sealing, as an adhesive in post and core applications, indirect bonding of veneers, inlays, onlays and crowns.

Order Information

31514 OptiBond Solo Plus Bottle Kit

Contains: 2 Bottles OptiBond Solo Plus (5 ml each), 1 Syringe Gel Etchant (3 g), 20 Gel Etchant Syringe Tips, 25 Four-Cavity Mixing Wells, 100 Applicator Tips

29682 OptiBond Solo Plus Unidose® Kit

Contains: 100 OptiBond Solo Plus Unidose Pouches (0.10 ml each – 10 ml total), 1 Syringe Gel Etchant (3 g), 20 Gel Etchant Syringe Tips, 100 Applicator Tips

31736 OptiBond Solo Plus Dual Cure Activator Bottle Kit

Contains: 1 Bottle OptiBond Solo Plus (5 ml), 1 Bottle OptiBond Solo Plus Dual Cure Activator (5 ml), 1 Syringe Gel Etchant (3 g), 20 Gel Etchant Syringe Tips, 25 Mixing Wells, 100 Applicator Tips

Refills

| 31513 | OptiBond | Solo | Plus | 15 ml hottlel |
|-------|----------|------|-------|-----------------|
| 31313 | Оривона | 2010 | 1 103 | (3 IIII boille) |

29669 OptiBond Solo Plus Unidose (pkg of 100 – 10 ml total)

31304 OptiBond Solo Plus Dual Cure Activator (5 ml bottle)

To learn more about our broad portfolio of products, visit KerrDental.com or contact your local Kerr representative.

OptiBond Family Overview

| | OPTIBOND FL | OPTIBOND SOLO PLUS | OPTIBOND ALL·IN·ONE |
|-------------------------|-------------|--------------------|------------------------|
| DIRECT APPLICATIONS | | | |
| INDIRECT APPLICATIONS | | * | ** |
| BOTTLE DELIVERY | | | |
| UNIDOSE DELIVERY | • | | |
| % FILL BY WEIGHT | 48 | 15 | 7 |
| Indirect Indications | | | |
| VENEERS | | | |
| PORCELAIN INLAYS/ONLAYS | | | |
| COMPOSITE INLAYS/ONLAYS | | | |
| CROWNS | | | |
| BRIDGES | | | • |
| CORE BUILDUP | • | * | • |
| Direct Indications | | | |
| COMPOSITE | | | |
| AMALGAM | | | |
| PORCELAIN REPAIR | • | • | |
| Bond Strengths (MPa) | | | |
| DIRECT TO DENTIN | 30 | 31 | 36 |
| DIRECT TO ENAMEL | 34 | 34 | 25 |
| INDIRECT TO DENTIN | _ | 29 | 30 |
| INDIRECT TO ENAMEL | _ | 26 | 32 |

^{*}With OptiBond Solo Plus Dual Cure Activator for dark cure applications.

Related Products

Whatever your preferred technique, total-etch or self-etch, Kerr has a solution for you— OptiBond Solo Plus and OptiBond FL® for total-etch or OptiBond® All•In•One for self-etch. For more information, visit KerrDental.com/bonding.









Data on file.

^{*} Indirect restorations with limited or no light accessibility—use OptiBond All-In•One with Maxcem™ self-adhering cement.